


Orbscan


Topographies import



Orbscan

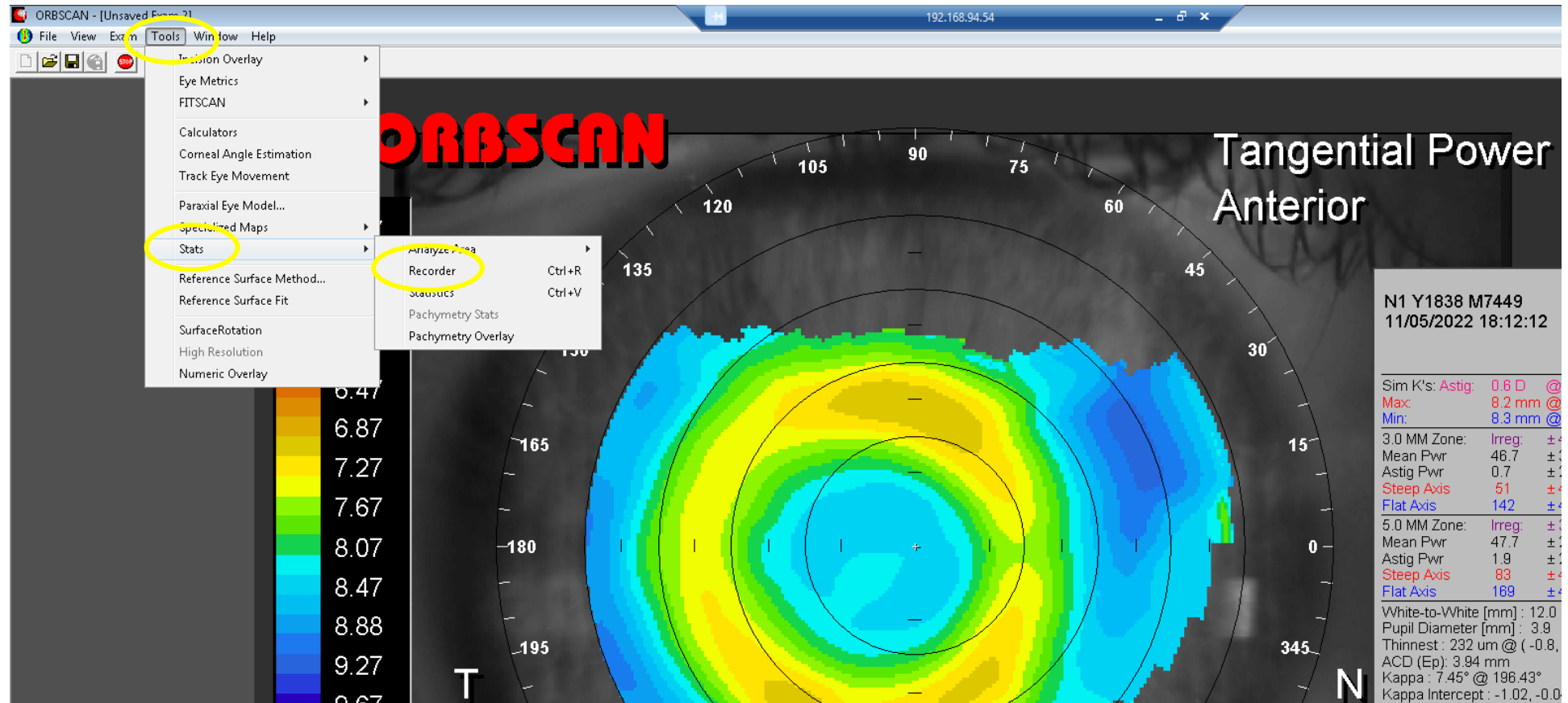

Orbscan

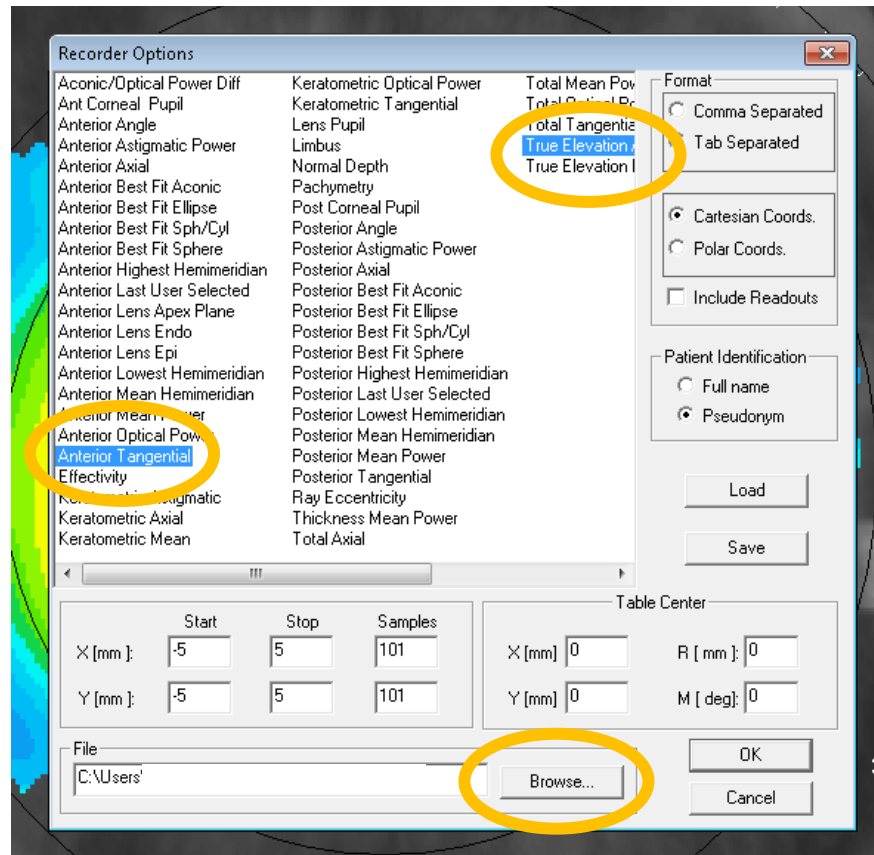

Save the exportations into the chosen folder and click on « OK »

Orbscan

Keratometry

Display unit
☒ Millimeters ☐ Diopters

Topographer model used
-

K flat	Exc flat	Axis (f)
	0,5	

K tight	Exc tight	DIV at 30° (Øcorneal)
	0,5	

K medium : 0.00

Topography





Save & import a topography for the Right Eye

Select a topographer*

Orbscan

ⓘ Attention, in the case of a pathological cornea:
Do not use the data cleaning function.



Clean topography



Import pivot format

For the Orbscan, one or two *.TXT files are expected

Select the file(s)

or drag and drop it here

Annuler

Import 

- In the patient's C&F folder, click on the symbol : 
- Drag the previously saved files into the frame that is displayed and click on

Importer 